


NO-DEVICE POSITION-OF-WEAR • MEASUREMENT GUIDE

1

Enable No Device and Position of Wear Modes

No Device Position of Wear mode allows you to capture all measurements without a device, including Vertex, Panto and Wrap.

Turn on the No Device and Position of Wear switches.

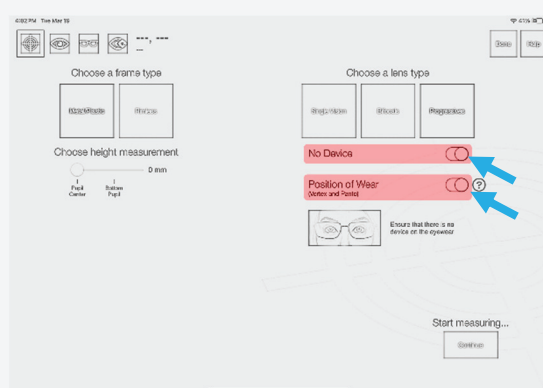

2

Prepare the Customer for the Picture Taking Process

Position the customer in their natural posture using standard opticianry techniques.

To ensure a smooth picture taking process, inform the customer of the upcoming picture taking steps.

Advise the customer to:


- Expect an iPad camera flash, followed by,
- An iPad down movement, followed by,
- A second iPad camera flash.

Ask them to follow the iPad with their eyes while holding their head still during the down movement and wait for the 2nd flash before changing their posture.



3

Capture the Measurement Images

- 1 • Ensure that the customer maintains their natural posture throughout the image capture process.
- 2 • Hold the tablet vertically in front of the customer's face while ensuring that the customer maintains their posture.
• Check the iPad-Tilt-Indicator to ensure the tilt is within the green zone.
- 3 • Ensure that the chin outline matches the outline of the customer's chin.
• Adjust distance between the iPad and the customer as well as the position of the iPad as necessary.
- 4 • Instruct the customer to look into the camera lens **while maintaining their head posture**.
• Tap the "Take Picture" button to capture the first image.

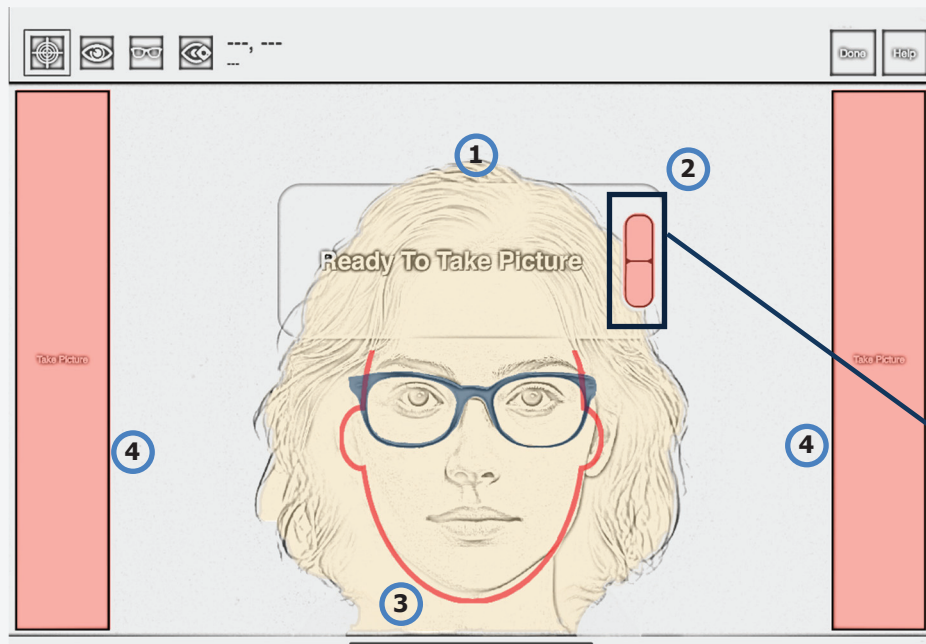

- 5 • Shift the iPad downwards until the new chin outline matches the outline of the customer's chin, while ensuring that the customer maintains their head posture.
• Repeat steps 1 to 4 to capture the second measurement image.

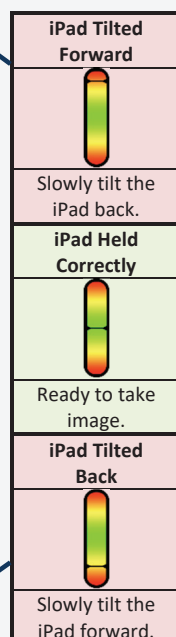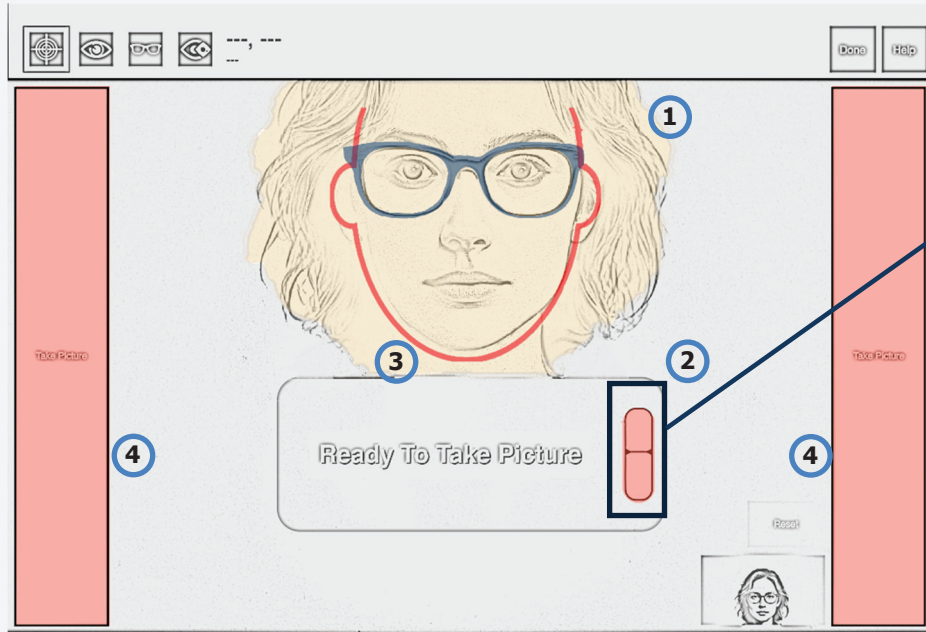


4

Verify the Marker Positions

Follow the on-screen steps and verify the position of the markers.


- 1 • Verify that the bridge marker is located at the top-center of the bridge in both images.
- 2 • Verify that the pupil markers are centered on the pupil reflection.
- 3 • Verify that the horizontal and vertical line markers are placed where the lens meets the frame in their respective images.

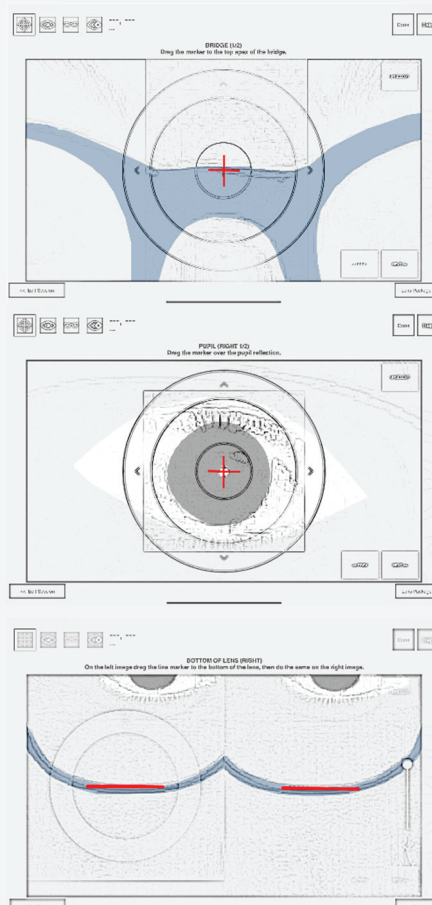

5

Measure the Frame Wrap

- 1 • Tap the **Tap to Take Wrap Measurement** text to begin the Wrap measurement session.

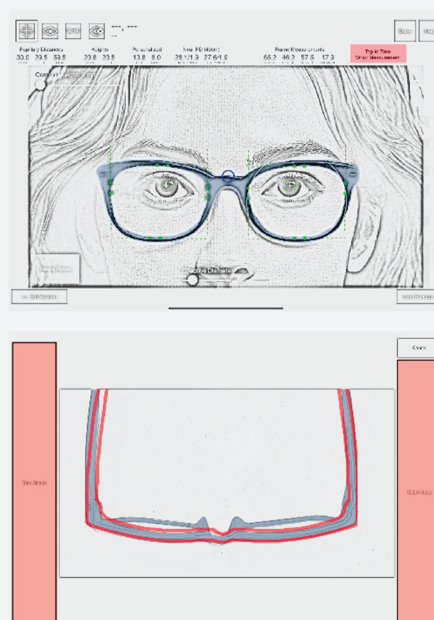

- 2 • Place the frame on top of an empty surface.
• Hold the iPad horizontally over the frame and adjust the distance so that the outline matches the outline of the frame.
• Capture the image when ready.

- 3
 - 1 • Place the Blue points at the outer edges of each lens.
 - 2 • Adjust the Wrap Angle slider so that the white line extends towards the inner edge of each lens.
 - 3 • Tap Unselect to complete the Wrap measurement.

