

OptikamPad DVI RxWizard Integration

OptikamPad integrates with DVI RxWizard. All measurements captured by the OptikamPad app can be transferred electronically to RxWizard, avoiding any manual entry of measurement results.

The Optikam Cloud portal at <https://cloud.optikam.com> allows the DVI Rx Wizard application to be launched with measurement fields pre-populated when the DVI Launcher is installed.

The integration requires the following:

- A Windows PC with DVI RxWizard and Optikam DVI Launcher application
- An active Optikam Cloud subscription

Steps to Configure:

- 1) Download and Install the Optikam DVI Launcher application from www.optikam.com/tools/dvilauncher.msi

Note: On older versions of Windows (prior to Windows 7 SP1) it may be required to install Microsoft .NET 4.0 (<http://www.microsoft.com/en-us/download/details.aspx?id=17851>)

Steps to perform measurement export:

- 1) Launch the Portal website: <https://cloud.optikam.com> from any browser and log into your Optikam cloud account.




Login

Username:

Password:

Login


- 2) In the Sessions tab select a saved session to export.

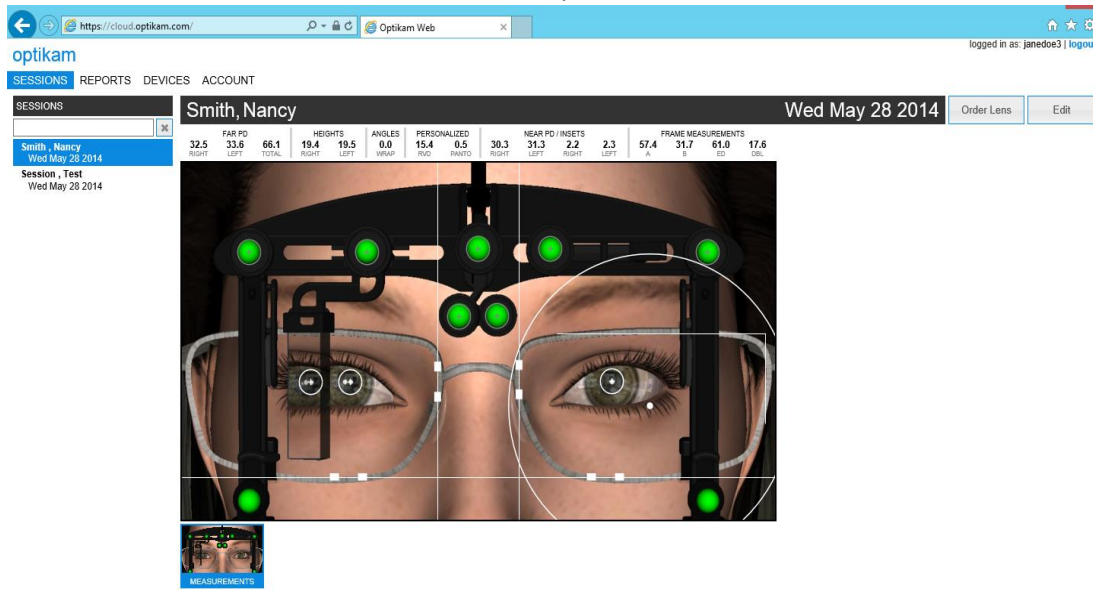

- 3) Click the “**Order Lens**” button to export the measurements to RxWizard.

The RxWizard application should now appear with the measurement fields populated.

If the RxWizard application was already running, the user will be presented the following popup instructing the user to close the RxWizard application and try again.

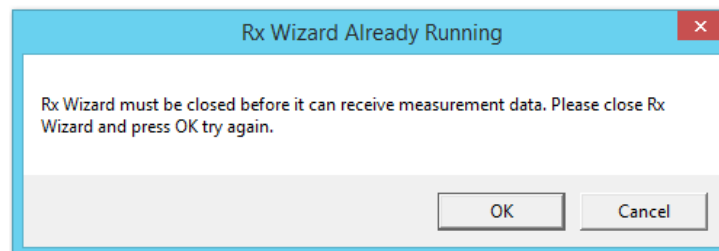

Note: If RxWizard is not installed in the default location (c:\dviremo) the following pop-up will prompt the user to locate the RxWizard installation directory.

